




# Baja Diffusive Air Blankets

## Warmer Patients, Better Outcomes

Baja Diffusive Air blankets are the warmest blankets ever in veterinary critical care. The patented 'Blue-side' material from AHS gently and uniformly delivers warm air to the patient without the disruptive flow from pin holes. If you're concerned about preventing hypothermia on longer cases, Baja's single-patient use Diffusive series has you covered pre-, intra- and post-op.


### Soft Tissue Blankets

- Convective Air, single-patient use
- Position over body, outside and under sterile field/draping
- Fenestrated to 'operate through'
- Available with, without tape around fenestration
- 4 sizes ensure most coverage on smallest, largest patients

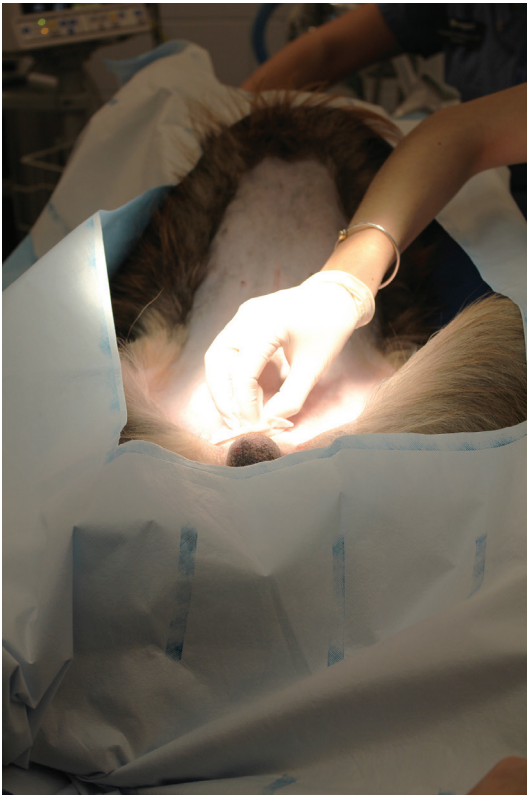

### Extremity Blankets

- Convective Air, 3 sizes
- Operative limb extends through fenestration
- Applied pre-prep; inflated in accordance with sterile technique
- Most patient coverage ever on difficult to drape extremity
- Outstanding normothermia maintenance on long cases

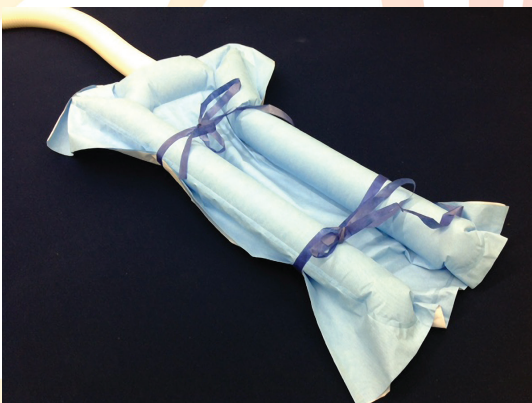

### Tube Blankets

- Underbody Blanket
- Most versatile Baja blanket
- Perfect for pre-op, intra-op or post-op warming
- Warms inwardly from perimeter U-shaped tube
- Non-woven ties – more fully wrap the blanket around the patient
- Two sizes to fit companion animals

*Baja blankets are designed for use with the Baja Desert Breeze blower*

